

Self-Collect Samples for HPV Testing

Available
at
Clinical
Labs

Women and people with a cervix from 25 to 74 years of age are encouraged to have a Cervical Screening Test every five years to check for the Human Papillomavirus (HPV). For Cervical Screening, patients have two options for HPV testing:

1. Clinician-Collected Samples

2. Self-Collected Samples

Self-collect HPV tests are available in-house and reported on at Clinical Labs for vaginal collection at home or within the medical clinic. Patients aged 25 to 74 are eligible for one bulk-billed HPV test every five years. Non-eligible tests will be privately billed.

Clinician-Collected versus Self-Collected samples for Cervical Screening tests:

	Clinician-collected cervical sample	Self-collected vaginal sample at home	Self-collected vaginal sample at clinic
Is co-testing possible?	Yes	No	No
Indicated for <ul style="list-style-type: none"> Those who are eligible and due or overdue for cervical screening, including during pregnancy Other points in the pathway where only a HPV test is required, including Test of Cure after HSIL treatment People who have had a total hysterectomy with history of HSIL 	Yes	Yes	Yes
<ul style="list-style-type: none"> Those undergoing surveillance after completely excised AIS 	Yes	No	No
<ul style="list-style-type: none"> People who have signs and symptoms suggestive of cervical cancer People who have postcoital or intermenstrual bleeding, post-menopausal bleeding, or unexplained persistent vaginal discharge People who were exposed to DES in utero 	Yes	No	No
<ul style="list-style-type: none"> STI Screening 	Yes	No	Yes



Patients should return the self-collect swab to their medical practice or a Clinical Labs collection centre with priority. Clinicians must forward self-collect samples to Clinical Labs within 21 days of collection – this ensures the sample is valid for testing. The Collection Date must be noted on the swab label.

Items required for self-collect at home:

- Clinical Labs request form signed by doctor (ensure self-collect HPV test at home is noted).
- Copan FLOQSwab® for vaginal self-collection (552C) (Red Top dry swab) – available from Clinical Labs.
- Only Copan FLOQSwab 552C is validated for use.




Items required for self-collect at clinic:

- Clinical Labs request form signed by doctor (ensure self-collect HPV test at clinic is noted).
- Copan FLOQSwab® for vaginal self-collection (552C.80) (Red Top dry swab) – available from Clinical Labs.
- ThinPrep® vial – available from Clinical Labs.
- Only Copan FLOQSwab 552C.80 is validated for use.



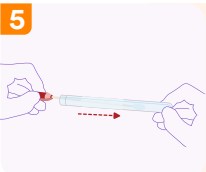
Steps for self-collecting at home:

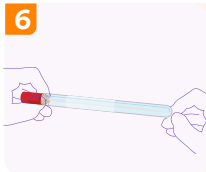
- 

Before you collect your sample, ensure you have clean hands. To avoid contamination, do not remove the swab from the container until you are ready to collect your sample and do not use any lubricant.
- 

Once you are comfortable, remove the red cap and pull the swab out of the container.
- 

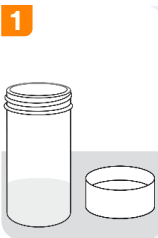
There is a red line halfway along the swab. Insert the front part of the swab inside your vagina aiming for this red line.
- 

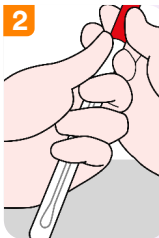
Once the swab is correctly inside your vagina, rotate it gently in a circular motion for 20-30 seconds.
- 


Remove the swab, place it back in the container and gently secure the cap.
- 

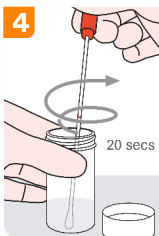
Return the self-collect swab to the medical practice or Clinical Labs collection centre with the request form with priority. Please ensure the Collection Date is labelled on the container.


Steps once the patient hands you the swab post-collection at medical practice:

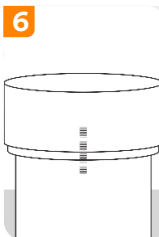
- 

Carefully uncup the ThinPrep® vial containing medium and place on a stable, flat surface.
- 

Slowly remove the FLOQSwab from the tube, taking care not to knock the swab against the walls of the tube.
- 

Fully immerse the FLOQSwab tip into the medium so it is touching the bottom of the ThinPrep® vial.
- 

Swirl the swab around the inner wall of the vial for 20 seconds, ensuring the swab remains immersed in the medium.
- 

Carefully draw the swab up the inner wall of the vial, draining any fluid from the swab. Discard the swab safely.
- 

Re-cap the ThinPrep® vial and tighten to prevent leakage. Store upright.
- Label ThinPrep® vial with the patient's full name and date of birth and complete the Clinical Labs patient request form. Please ensure **self-collect HPV test (at medical practice)** is noted on the request form.
- Send the ThinPrep® vial to the Clinical Labs laboratory, as per normal process.

About the author:



Dr Marissa Muller

MBChB (UL) MMed (Path) Anatomical Path.(UP)

Lab: Subiaco

Speciality: Anatomical Pathology

Areas of Interest: Hematolymphoid, Gastrointestinal, Soft tissue, Head and neck, Lung and Pediatric pathology

Email: marissa.muller@clinicallabs.com.au



For more information on cervical screening, visit clinicallabs.com.au/cervicalscreeningprogram